


TO ADJUST RANGE OF MOTION (ROM):

- Step 1:** Press the circle with the arrow pointing upward 1 x to make the number on the number on the right side of the window flash on and off.
- Step 2:** Only when the number is flashing can ROM be increased by pressing the triangle with the (+) SIGN.
- Step 3:** Looking at the window on the remote, you should see the number to the right increasing as you press the (+) sign.
- Step 4:** Press the START button once ROM has been set to begin treatment session.
- Step 5:** Press the STOP button when session is complete.

TO ADD ROTATION:

- Step 1:** Begin by pressing the ADD/ABD WITH ROTATION button at the top right of the remote until you see 'ROT' flashing in the window.
- Step 2:** Repeat steps 1 - 5 above.

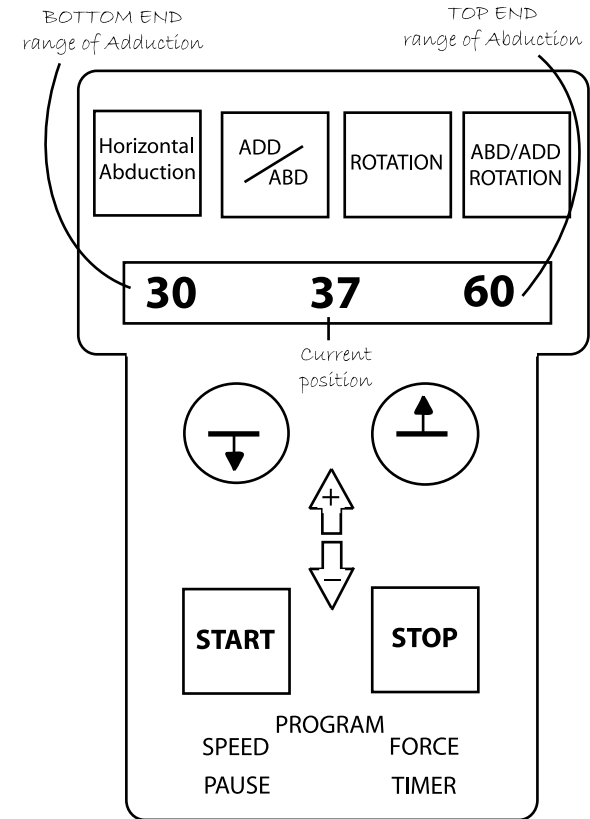

PLEASE CALL **1-800-339-9295**

EXT. 301 (LA) 302 (SBN/RIV) 303 (OC) 304 (VENTURA) 305 (SD)
FOR QUESTIONS, PROBLEMS, AND PICK UP OF EQUIPMENT